

## Downloading and installing the free open-source database postgresSQL for deltaMasses

Time needed : approx. 10 minutes.

You might want to contact your system administrator for assistance.


Open <http://www.postgresql.org/download/windows/> and click „download“

which takes you to


<http://www.enterprisedb.com/products-services-training/pgdownload#windows>

Assuming you have a Windows 64 bit machine (which you must have for deltaMasses), download version 9.2.2 (or better) for Win X86-64. Size 50.9 Mbyte.

You will get the file postgresql-9.2.2-1.windows-x64.exe. Start installing.

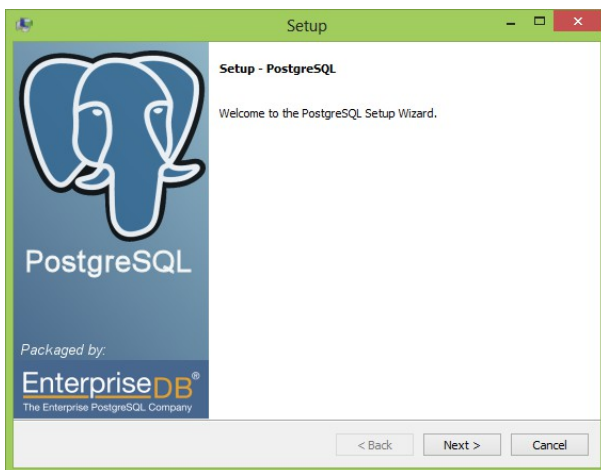


->Next

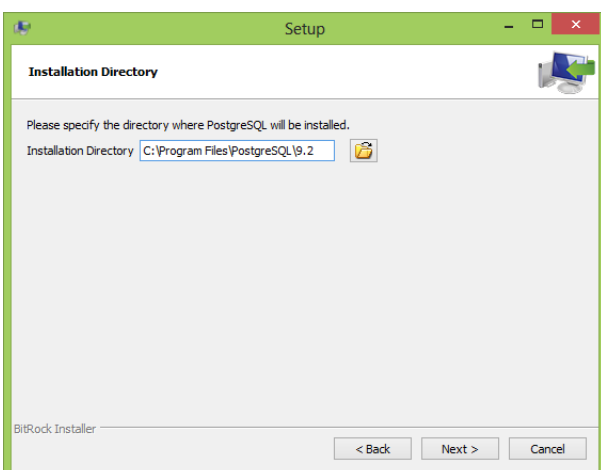


-> your choice. The default suggested by the program is typically OK. If you have a SSD , you might want to put it there for speed reasons.

-> Next

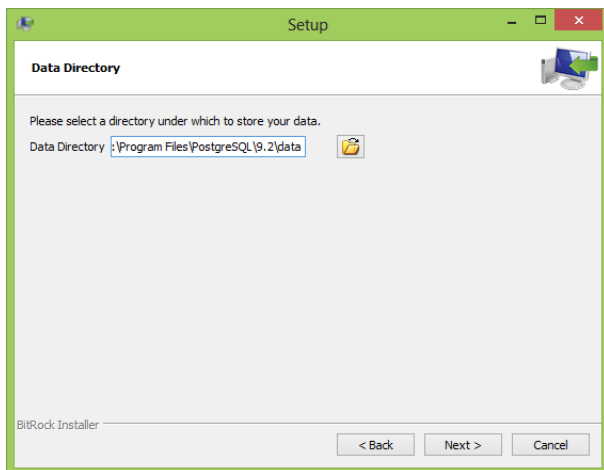


Where the data is stored. (For those who have an SSD drive: put it there, runs faster). Most users: leave it as suggested.

-> Next

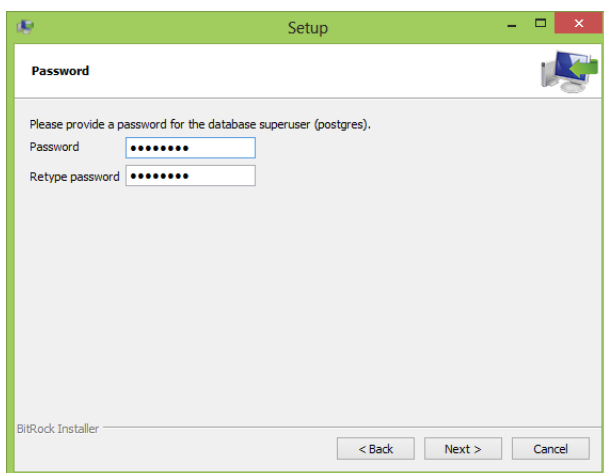


Enter your database password of choice ( PLEASE REMEMBER THIS PASSWORD ) The username is postgres. We refer to this password below by DB\_PASSWORD . Do not use a password you use elsewhere. Our default password ist 4.3.jjMM

-> Next

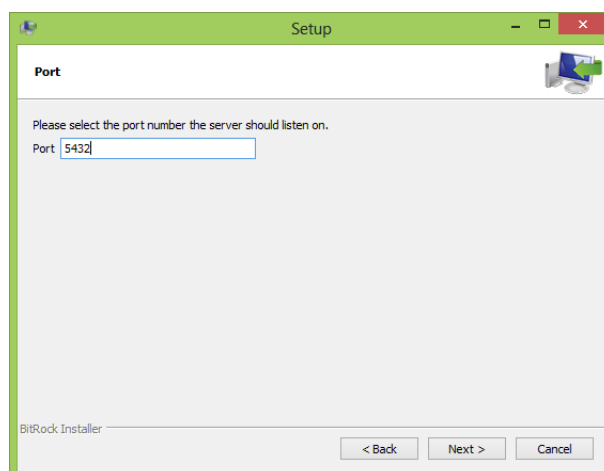


Port configuration. You will typically leave it as suggested. PLEASE REMEMBER YOUR PORT NUMBER. In most cases you will get the 5432.

-> Next

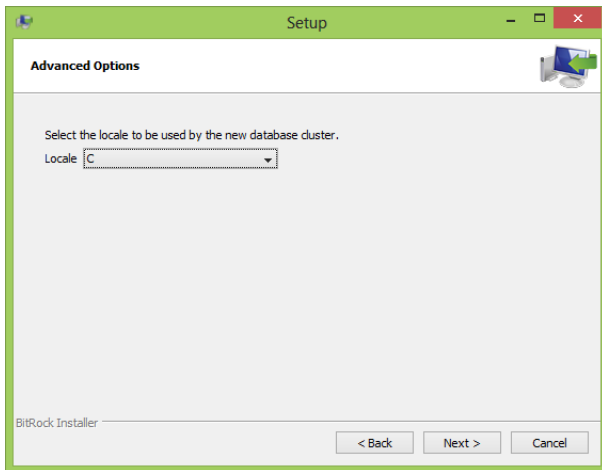


Set the Locale of the database to „C“. -> Next

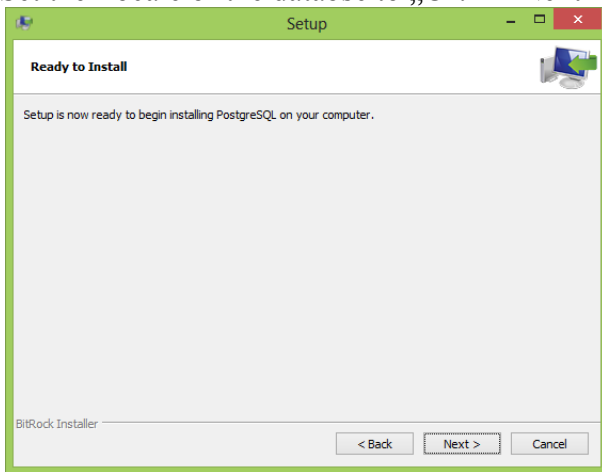

-> Next

The installation will run a couple of minutes. On our system it takes about 3 minutes.

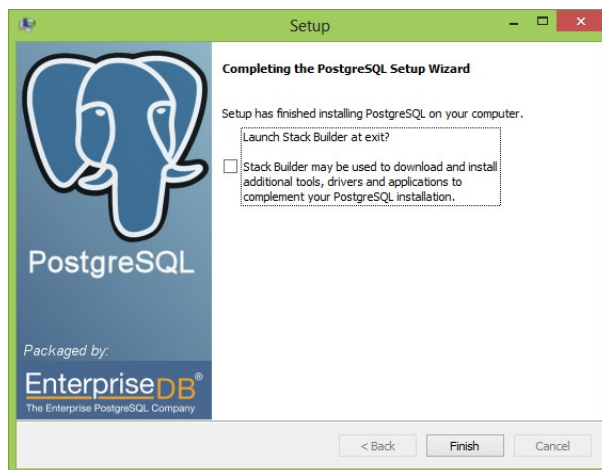

-> Finish.